




TAB ADHESIVE ELECTRODES

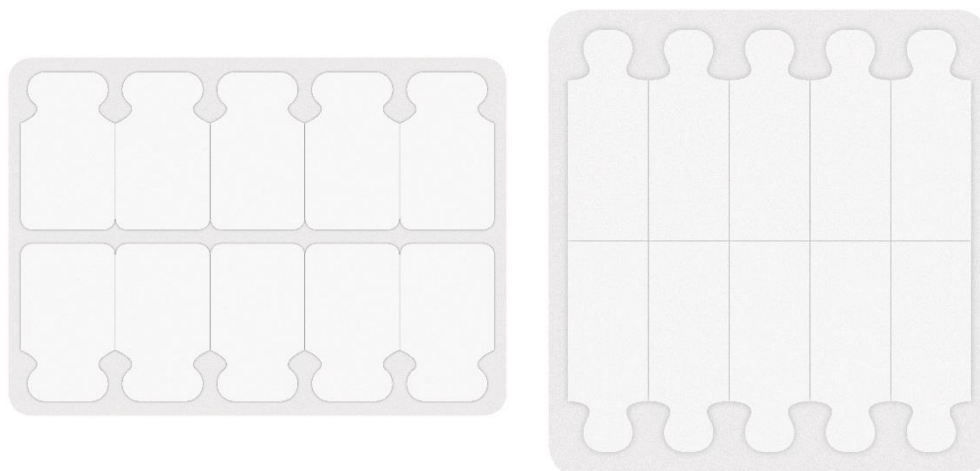

Description

Spes Medica adhesive electrodes have a special shape and a specific size studied to allow a perfect registration even for the smaller muscles.

The solid gel is developed especially for this application and gives to the electrode a stable adhesion on the skin. It also can be repositioned several times.

The electrode is made by biocompatible soft materials that give the best comfort to the patient.

- High quality gel for the best adhesiveness;
- AgAgCl electrode for the best conduction;
- Easy tab connection with all clip alligator connectors.

UMDNS 11440

Materials

Safety

Specifications

Environment

Electrode	AgAgCl
Material in contact with the skin	Hydrogel
Liner	PET
Packaging	AL/PE

Biocompatibility (UNI EN ISO 10993)	Yes
Latex	No
Phthalates	No

Generic code	TABSMXXXXXXX
Dimension	14x34 mm 20x34 mm
Connector	Tab

RoHS II	✓
REACH	✓